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INSTRUCTIONS

Use of the electric toothbrush

Hold the brushhead under running water for a few seconds and then apply some toothpaste.

Do not switch the toothbrush on until it is in your mouth. This will ensure that you avoid splattering the taps, mirror and your clothing. Press the upper part of the switch "1" to switch the toothbrush on (A).

When brushing, keep the brush parallel to the teeth (B). Guide the brush slowly over the teeth and the edge of the gums, brushing each tooth individually (C). First brush the outside faces of the teeth (D), then the inside faces (E) and finally the chewing surfaces (F).

Brush the teeth using the toothbrush for at least two minutes to ensure they are thoroughly and effectively cleaned.



For an even cleaning effect, change the direction of cleaning each time of use. Cleaning your teeth regularly, morning and evening, like this prevents build-up on the teeth and keeps teeth and gums healthy.

For the first few days of use the toothbrush may cause some bleeding of the gums. This only lasts until the gums get used to the oscillating movement of the brushhead. If this bleeding continues for more than two weeks, then you should consult a dentist.

After brushing the teeth switch off the toothbrush using bottom part of the switch "0" and only then take the brush out of your mouth. Then clean the brushhead and the handle thoroughly with the toothbrush switched off.

Shake off any superfluous water and dry the brush with a dry towel.

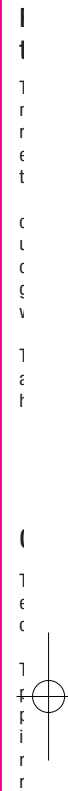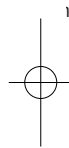Replacing the brushhead

To ensure that the toothbrush maintains its cleaning performance, we recommend that the brushhead is changed every three months and certainly by the time the bristles start to weaken (G). If you find you have to change your brushes more often than every three months, then you are using too much pressure when cleaning your teeth. Always let the brush glide smoothly over the teeth without using too much pressure.

To change the brush, pull out the used one and press the new brush firmly into the handle.

Changing the batteries

To obtain the best performance from the electric toothbrush we recommend the use of AA NiCd batteries.

To change the batteries, remove the bottom part of the handle (grip in the recesses and pull outwards) (H). Slide the new batteries into the handle and use firm pressure to replace the bottom part. You can use rechargeable batteries in your toothbrush.

Remove the batteries if the toothbrush is not to be used for a long time. Use 2 AA batteries.

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E
E

F
C
H

- ⓐ
- A. Press the switch to switch the toothbrush on.
 - B. When brushing, keep the brush parallel to the teeth.
 - C. Guide the brush from one tooth to the next over the space between the teeth.
 - D. Thoroughly clean the front and back faces of your back teeth.
 - E. Clean the inside faces of your front and canine teeth.
 - F. Then finally brush the chewing surfaces.
 - G. Change the brushheads every three months.
 - H. Change the batteries when required.

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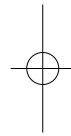
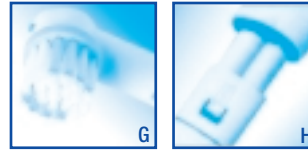
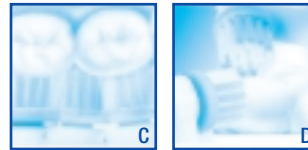
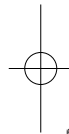
(GB)

ELECTRIC TOOTHBRUSH

- Ideal for travel
- batteries included

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