SILVER CREST®

www.lidl-service.com

SHIATSU NECK MASSAGER SSMN 2 B3

GB IE NI

SHIATSU NECK MASSAGER

Operating instructions

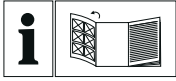
DE AT CH

SHIATSU-NACKENMASSAGE- HÖRNCHEN

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IAN 290496

IE NI

GB IE NI

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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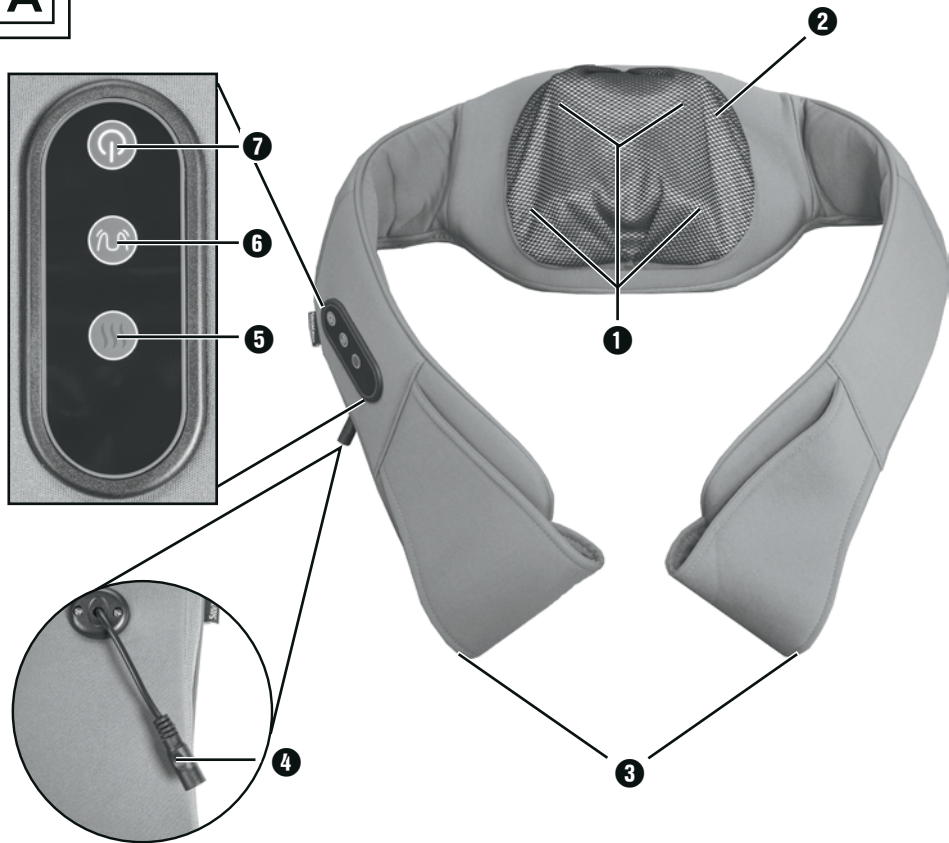
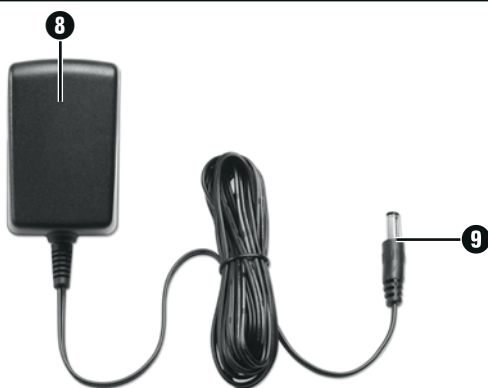
A**B**

Index

Introduction	2
Copyright	2
Intended use	2
Items supplied	3
Disposal of the packaging	3
Appliance description	4
Technical data	4
Safety instructions	5
What is Shiatsu massage?	7
Operation	7
Connection	7
Fitting	8
In use	8
After the massage	9
Cleaning	9
Cleaning the removable cover	9
Cleaning the appliance	10
Storage	10
Troubleshooting	11
Disposal of the appliance	11
Disposal of the mains adapter	11
Kompernass Handels GmbH warranty	12
Service	13
Importer	13

Introduction

Congratulations on the purchase of your new appliance.

You have clearly decided in favour of a quality product. These operating instructions are a part of this product. They contain important information in regard to safety, use and disposal. Before using the product, familiarise yourself with all of these operating and safety instructions. Use the product only as described and only for the specified areas of application. Retain these instructions for future reference. In addition, pass these documents on, together with the product, to any future owner.

Copyright

This documentation is copyright protected.

Any copying or reproduction of it, including as extracts, as well as the reproduction of images, also in an altered state, is only permitted with the written authorisation of the manufacturer.

Intended use

This appliance is intended for massaging the neck and shoulders of people.

This appliance is intended solely for personal use in the domestic sector, it is not to be used for medical or commercial applications and may not substitute for medical treatment. The appliance is designed for use in dry indoor rooms. It is not intended for use in the medical/therapeutic or commercial fields.

WARNING

Risks from unintended use!

Risks can be engendered if this appliance is used for purposes not intended for it and/or for other types of use.

- ▶ Use the appliance exclusively for its intended purposes.
- ▶ Pay heed to the procedures described in these operating instructions.

NOTICE

- ▶ Risks can be engendered if this appliance is used for purposes not intended for it and/or for other types of use. Use the appliance exclusively for its intended purposes. Observe the procedures described in these operating instructions. Claims of any kind for damage resulting from unintended use, incompetent repairs, unauthorised modification or the use of unauthorised spare parts will not be accepted. The operator alone bears liability.

Items supplied

The appliance is delivered with the following components as standard:

- Shiatsu Neck Massager
 - Mains adapter
 - Operating Instructions
- 1) Remove all parts of the appliance and the operating instructions from the carton.
 - 2) Remove all packaging material.
 - 3) Clean all parts of the appliance as described in the section "Cleaning".

NOTICE

- ▶ Check the contents to ensure that everything is present and for signs of visible damage.
- ▶ If the delivery is not complete or something is damaged due to defective packaging or through transportation contact the Service Hotline (see chapter **Service**).

Disposal of the packaging

The packaging protects the appliance from transport damage. The packaging materials are selected from the point of view of their environmental friendliness and disposal technology and are therefore recyclable.

The recirculation of packaging into the material circuit saves on raw material and reduces the amount of waste generated. Dispose of packaging material that is no longer needed as per the regionally established regulations.

NOTICE

- ▶ If possible, preserve the appliance's original packaging during the warranty period so that, in the event of a warranty claim, you can package the appliance appropriately for its return.

Appliance description

- ❶ Massage heads
- ❷ Removable cover
- ❸ Hand loops
- ❹ Socket
- ❺ "Warmth" button 
- ❻ "Massage" button 
- ❼ "On/Off" button 
- ❽ Mains adapter
- ❾ Connection plug

Technical data

Mains adapter ❽	
Type	YJS010A-1201000B
Input voltage	100 - 240 V ~, 50/60 Hz
Current consumption	350 mA
Output voltage	12 V  (Direct current)
Output current	1000 mA
Polarity	
Safety transformer	
Switched-mode power supply	
Nominal ambient temperature	T _a =35°C
Protection type	IP20 Protection against solid objects with a diameter of more than 12.5 mm
Appliance	
Input voltage	12 V 
Current consumption	1 A
Polarity	

Safety instructions

⚠ RISK OF ELECTRIC SHOCK

- ▶ Operate the appliance only with the supplied power adapter YJS010A-1201000B.
- ▶ Connect the appliance only to correctly installed and earthed mains power sockets. Ensure that the rating of the local power supply tallies completely with the details given on the rating plate of the appliance.
- ▶ Ensure that the power cable cannot become wet or damp while the appliance is in operation. Place it such that it cannot be crushed or damaged.
- ▶ Keep the power cable away from hot surfaces.
- ▶ Do not perform any repairs on the appliance. All repairs must be performed by our Customer Services or by qualified personnel.
- ▶ Always remove the plug from the mains power socket before you clean the appliance.
- ▶ Do not use the appliance in damp areas, such as in the bathroom.
- ▶ To avoid risks, arrange for defective plugs and/or cables to be replaced at once by qualified technicians or our Customer Service Department.
- ▶ Do not insert needles or other sharp objects into the appliance.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not use an external timing switch or a separate remote control system to operate the appliance.
- ▶ This appliance may be used by children aged 8 over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use of the appliance safely and are aware of the potential risks. Children must not play with the appliance. Cleaning and user maintenance tasks must not be carried out by children unless they are supervised.
- ▶ Do not massage infants, young children, helpless persons or animals with the appliance.
- ▶ NEVER use the appliance if you have injuries or pain in the neck or cervical spine.
- ▶ Stop treatment with the appliance IMMEDIATELY if you feel it to be uncomfortable or painful.
- ▶ During use the appliance has a hot surface. Individuals who are insensitive to heat must take care when using the appliance.
- ▶ Do not use the appliance on swollen, burnt, inflamed, diseased or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the appliance.
- ▶ Ensure that your fingers or other body parts do not come between the rotating massage heads. Risk of crush injuries!

⚠ WARNING! RISK OF INJURY!

- ▶ The appliance may not be used in the face (e.g. eyes), on the throat, feet or other sensitive body parts.
- ▶ During use do not cover the appliance (with a pillow, blanket, etc.). NEVER use the appliance in the vicinity of petrol or other inflammable materials.

WARNING - PROPERTY DAMAGE!

Use the appliance only indoors.

- ▶ Operate the appliance only with the supplied power adapter.

What is Shiatsu massage?

Shiatsu (shi = finger, atsu = pressure) is a form of massage developed in Japan, which itself emerged from traditional Chinese medicine. Shiatsu massage consists of a relaxing and powerful massage with fingers, thumbs or palms. This massager mimics the movements of shiatsu massage with the help of rotating massage heads ❶.

Operation

Connection

- 1) Connect the connection plug ❹ on the mains adapter ❸ with the socket ❷ on the appliance.
- 2) Connect the mains adapter ❸ into a mains power socket.

Fitting

- 1) Place the appliance around your neck so that the massage heads **1** lie on the sides of the neck. You can use the hand loops **3** to align it.
The control panel with the buttons must be on your right.

In use

Preferably, lie back for the massage and relax. Provide for a relaxing atmosphere. For example, dim the lights a little.

NOTICE

- Ensure that you experience the massage as comfortable and relaxing. If needs be, vary the pressure and/or your position. Stop the massage immediately if it is painful or uncomfortable.

NOTICE

- Never use the appliance for more than 15 minutes. A longer massage may lead to tension caused by overstimulation of the muscles.

You can select between 2 functions:

- Shiatsu massage with the massage heads **1** (fast and slow)
- Warmth (2 of the massage heads **1** illuminate and warm up)

You can combine various massage functions with each other or use them individually:

- Press the "On / Off" button  **7**. The massage heads **1** begin with the slow massage.
- Press the button "Massage"  **6** to increase the massage speed. Press the button "Massage"  **6** once again to stop the massage heads **1**.
- Press the "Warmth" button  **5**, to switch the warmth on. The "Warmth" button  **5** and the massage heads **1** light up. Press the "Warmth" button  **5** again to switch the warmth function off. The light on the massage heads **1** as well as the button "Warmth" button  **5** goes out.

NOTICE

- Once you have switched the warmth function on and the massage heads **1** are rotating and you then stop the massage heads **1** by pressing "Massage"  **6**, the warmth function also turns off at the same time. If you want to use the warmth function again, press the "Warmth" button  **5** to reactivate the warmth function.

NOTICE

► The appliance switches off automatically after 15 minutes.

- If you want to terminate the massage before the end of 15 minutes and wish to switch the appliance off, press the "On/Off" button  **7**.

After the massage

- 1) Disconnect the connection plug **9** from the socket **4**.
- 2) Remove the appliance.
- 3) Disconnect the mains adapter **8** from the mains power socket.

Cleaning

Cleaning the removable cover

- Switch off the appliance and then pull the mains adapter **8** out of the socket.
- Completely open the zip on the removable cover **2** to take the removable cover **2** off the appliance.
- Note the following washing instructions when you wash the removable cover **2** in the washing machine:

	max. 30°C in the delicates cycle
	do not bleach
	do not tumble dry
	do not iron
	do not dry clean

The cover material is 70% nylon and 30% polyester.

Cleaning the appliance

RISK OF ELECTRIC SHOCK

- ▶ Before you clean the appliance, always disconnect the mains adapter  from the mains power socket.
- ▶ UNDER NO CIRCUMSTANCES should you immerse the appliance in water or other liquids during cleaning!

WARNING - PROPERTY DAMAGE!

- ▶ Do not use any aggressive cleaning liquids or solvents. They could damage the surface.

NOTICE

- ▶ The zipper on the appliance must not be opened, because it is fitted only for production-technical reasons.

- 1) Wipe the appliance with a damp cloth. If necessary, add a little mild detergent onto the cloth.
- 2) Wipe the mains adapter  off with a damp cloth.
- 3) Let all parts dry thoroughly before re-use.

Storage

- Clean the appliance as described in the chapter "Cleaning".
- Store the appliance at a clean and dry location.

Troubleshooting

Defect	Cause	Remedy
The appliance does not work.	The mains adapter ⑧ is not connected into a mains power socket.	Connect the mains adapter ⑧ with a mains power socket.
	The appliance is not switched on.	Switch the appliance on.
	The appliance is defective.	Contact Customer Services.
	The connection plug ⑨ is not connected to the socket ④ on the appliance.	Connect the connection plug ⑨ with the socket ④.

Disposal of the appliance

Do not dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance through an approved disposal centre or at your community waste facility. Pay heed to the currently applicable regulations. In case of doubt, please contact your waste disposal centre.

Disposal of the mains adapter

Do not dispose of the mains adapter with the normal household waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the mains adapter at an approved disposal centre or your community waste facility. Please observe applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on www.lidl-service.com.

Service

GB Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: kompennass@lidl.co.uk

IE Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: kompennass@lidl.ie

IAN 290496

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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