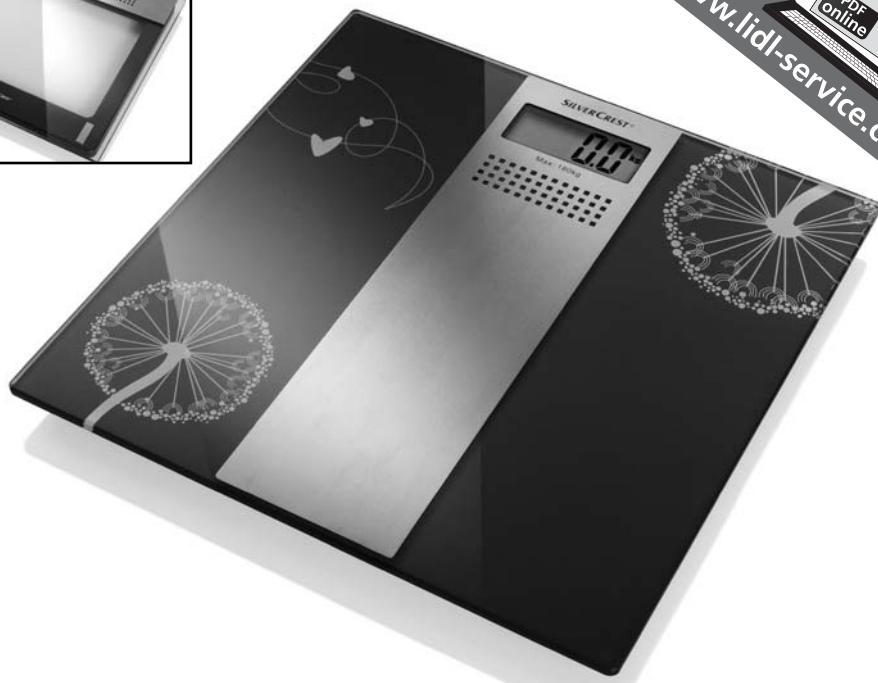
# SILVERCREST®

[www.lidl-service.com](http://www.lidl-service.com)

## TALKING SCALE SPWS 180 B2

(GB) (CY)

### TALKING SCALE

Operating instructions

(RO)

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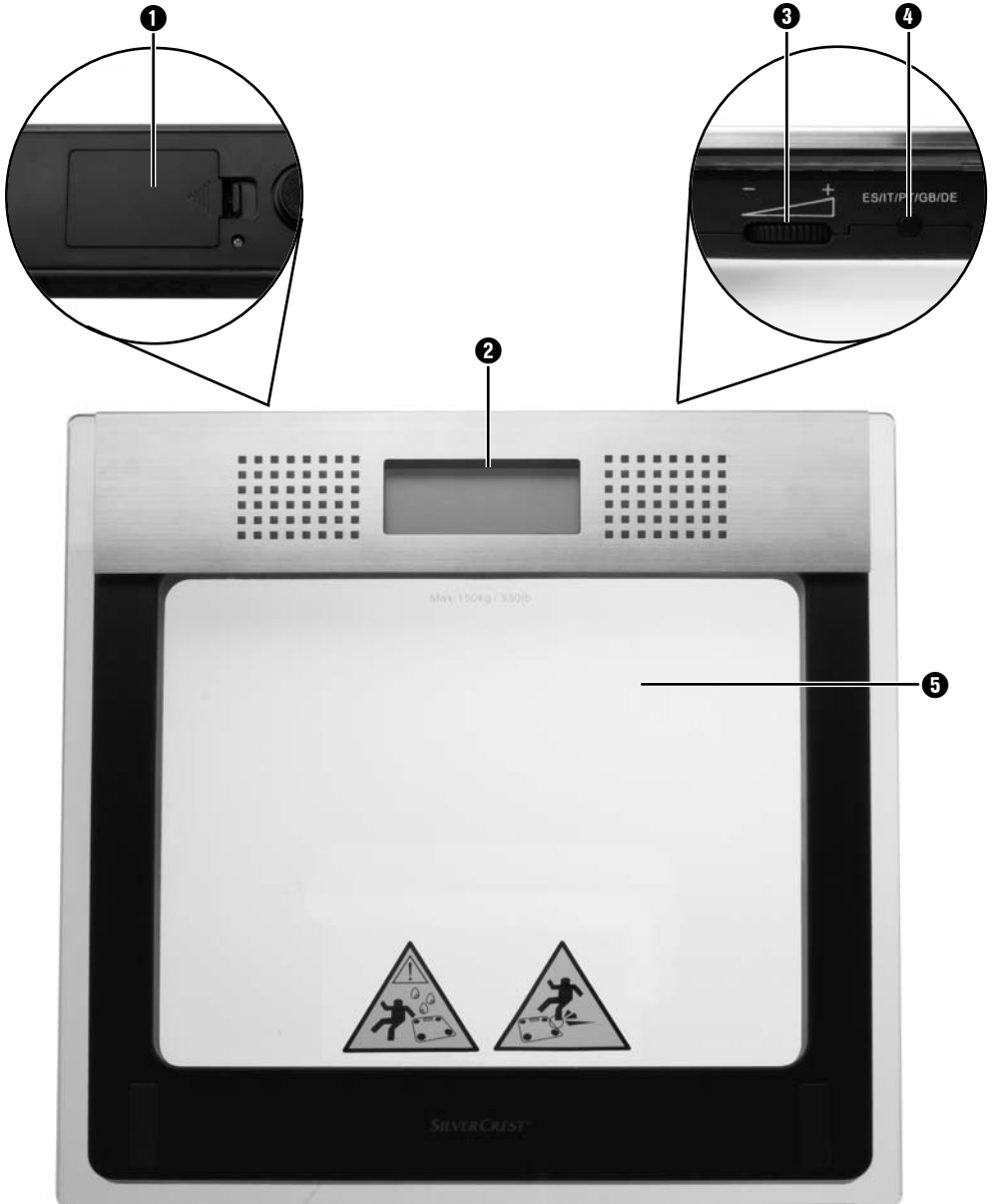
<b>Index</b>	<b>Page</b>
<b>Overview of the appliance</b>	<b>2</b>
<b>Introduction</b>	<b>4</b>
<b>Intended Use</b>	<b>4</b>
<b>Items supplied</b>	<b>4</b>
<b>Description of the appliance</b>	<b>4</b>
<b>Technical Data</b>	<b>4</b>
<b>Safety instructions</b>	<b>5</b>
<b>Using the scale</b>	<b>7</b>
<b>Warning indicators</b>	<b>8</b>
<b>Replacing the batteries</b>	<b>8</b>
<b>Usage and safe-keeping</b>	<b>8</b>
<b>Cleaning</b>	<b>8</b>
<b>Troubleshooting</b>	<b>9</b>
<b>Disposal</b>	<b>9</b>
<b>Warranty and Service</b>	<b>10</b>
<b>Importer</b>	<b>10</b>



Read the operating instructions carefully before using the device for the first time and preserve this booklet for later reference. Pass this manual on to whoever might acquire the device at a future date.

# Talking Scale

## Overview of the appliance

Model 1

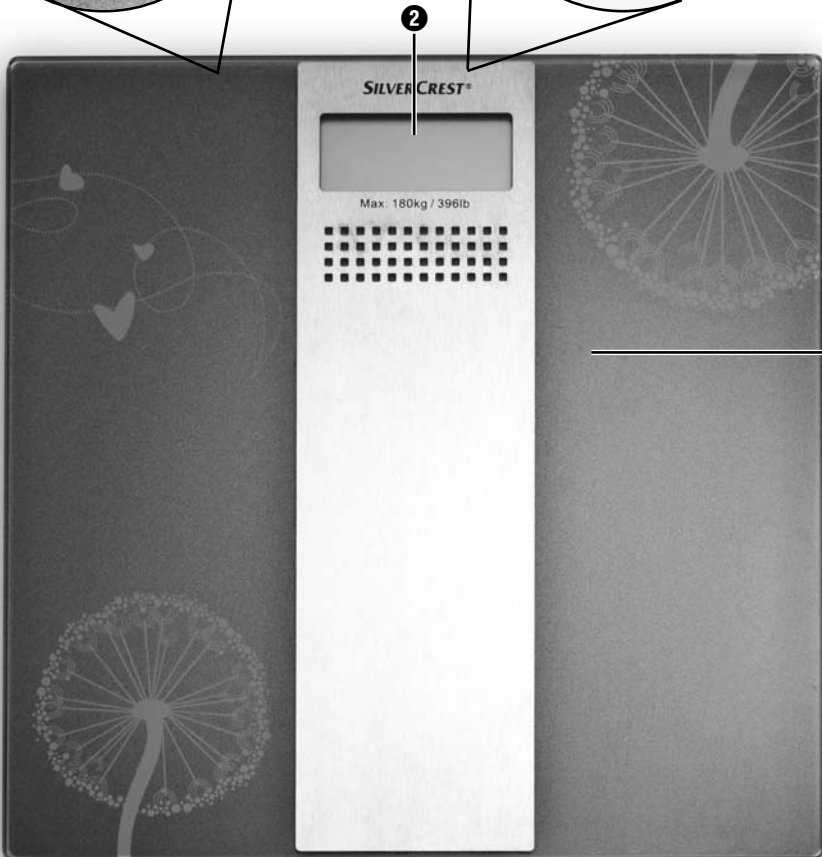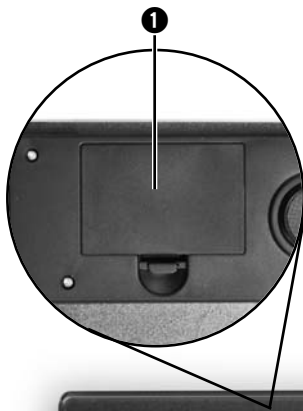Model 2

## Introduction

Congratulations on the purchase of your new appliance.

You have clearly decided in favour of a quality product. These operating instructions are a part of this product. They contain important information in regard to safety, use and disposal. Before using the product, familiarise yourself with all of these operating and safety instructions. Use the product only as described and only for the specified areas of application. In addition, pass these documents on, together with the product, to any future owner.

## Intended Use

This scale is intended solely

- for use as a bathroom scale
- in a domestic environment.

This appliance is not intended for commercial or industrial use.

## Items supplied

Talking Scale

3 x Batteries Type AAA/Micro/LR03

Operating instructions

Before using the appliance for the first time, check to ensure that the contents are complete and undamaged.

## Description of the appliance

- 1 Battery compartment
- 2 Display
- 3 Volume control
- 4 Language button
- 5 Safety glass

## Technical Data

Maximum weighing capacity: 180 kg / 396 lb

Tolerance range: +/- 3%

Minimum measurable

weight/Minimum Weight: 3 kg / 6.6 lb

Division: 0.1 kg / 0.2 lb

Operating temperature: 10°C - 40°C

Power supply : 4.5 V  =

Batteries: 3 x 1.5 V DC  
(Type AAA/Micro/LR03)

Display for weak battery

Overload indicator

Error indicator

Abbreviations for the selected display languages **2**:

Display	Language
GB	English
HR	Croatian
RO	Romanian
GR	Greek
BG	Bulgarian

## Safety instructions

---

### **⚠ *Danger!***

- Ensure that no liquids permeate into the scale, the display or the speaker. There is a risk of receiving an electric shock!
- This appliance may be used by children aged 8 over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks. Children must not play with the appliance. Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.

### **⚠ *Warning!***

- Do not place the scale in the immediate vicinity of sources of heat. The scale and the display may only be operated at temperatures between 10° and 40°C. Higher or lower temperatures can lead to malfunctions or to permanent damage.
- Never burden these bathroom scales with a weight of more than 180 kg / 396 lb! This could irreparably damage the scales.

### **⚠ *Information regarding the handling of batteries***

This appliance uses batteries. For your interaction with batteries, please take note of the following:

### ***Risk of explosion!***

Do not throw batteries into a fire. Do not recharge the batteries.

- Never open batteries, never solder or weld batteries.  
The risk of explosions and injuries exists!
- Regularly check the condition of the batteries. Leaking batteries can cause damage to the appliance.
- If you do not intend to use the appliance for an extended period, remove the batteries.
- In the event of the batteries leaking acids, put on a pair of protective gloves. Then clean the battery compartment and the battery contacts with a dry cloth.
- In the event of skin contact with leaked battery acid, immediately rinse the affected area with lots of clean water.
- Keep batteries away from children. Children can put batteries into their mouths and swallow them. If a battery is swallowed, medical assistance must be sought immediately.

**ⓘ Notice!**

Avoid electrostatic charging and discharging. These can cause erroneous displays and/or the sudden deactivation of the scales.

For damage caused by...

- opening of the scales housing
- repairs attempted by non-qualified personnel and
- improper usage of the scales

no liability/warranty will be accepted by the manufacturer!

## Using the scale

1. Open the battery compartment ❶. Insert the batteries as per the polarity diagram in the battery compartment ❶.
2. Close the battery compartment ❶. The battery compartment lid must audibly engage.
3. Place the scale on a **flat and stable surface**.  
Ensure that neither the floor at the location of the scale nor the upper surface of the scale are wet, otherwise there is a risk of slipping!

### ❶ **Notice!**

This scale contains sensitive electronic components and it should therefore NOT be located:

- close to sources of heat
- in environments with extreme air humidity.

#### **Attention:**

Never burden the scales one-sidedly. Always step onto the scales in the middle and with both feet. Otherwise, the scales could tip and become irreparably damaged. Additionally, there is the risk of falling over!

#### **Attention:**

Never use the scale if it is wet or moist. Never step onto the scale with wet or moist feet. There would be a serious risk of slipping and personal injury!

### ⚠ **Attention:**

NEVER jump onto the glass scale, this could irreparably damage it.

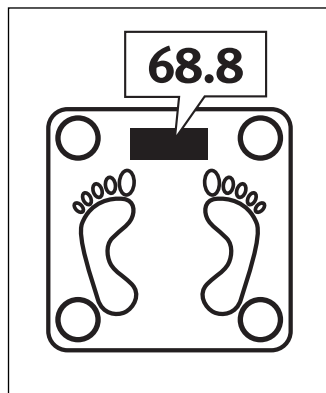
4. Stand briefly on the scales so that the device turns itself on.

5. Press the language button ❷ repeatedly until the desired language appears on the left of the display ❸.
6. Set the desired volume level using the volume control ❹. To deactivate the voice, turn the volume control ❹ towards the "0" until you hear a click.

### ❶ **Notice!**

The language and the sound volume must only be adjusted on the first use.

7. Stand with both feet on the scales. Stand as still as possible and without displacing your weight.
8. During the weighing process (about 5 seconds) the displayed value increases until your weight is determined.
9. When your weight is determined, the result blinks twice on the display ❸ and it speaks your determined weight.

After you have stepped down from the scale it sounds "Ready for operation" and you can then start another measurement.

If you simply wait, the scale will indicate by a double beep that it has switched itself off.

## Warning indicators

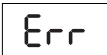
---

Should your weight exceed the maximum weighing capacity of 180 kg / 396 lb, "0 - Ld" is shown in the display ②. You should IMMEDIATELY step down from the scale, otherwise it could be damaged.

If the display ② shows "LO" and "Low Battery" sounds, you must replace the batteries.

If the display ② shows "Err" and "Error" sounds, the scale cannot measure. In this case, step down from the scale, wait a few seconds and then try again. Take note that you need to both stand still and not interrupt the weighing procedure.

## Replacing the batteries

---

1. Open the battery compartment ❶ and remove the batteries.
2. Insert 3 new batteries of the type Micro/AAA/LR03 as per the polarities specified in the battery compartment ❶.
3. Close the battery compartment ❶. The battery compartment lid must click in audibly.

## Usage and safe-keeping

---

- Avoid unnecessary battery discharge. Accordingly, do not place objects on the glass scale when it is not being used. Remove the batteries if you do not intend to use the scale for an extended period.
- Keep the scale at a clean and dry location.
- Avoid electrostatic charging.
- Your scale is a sensitive electronic device. As such, it can occasionally be disturbed by other high-frequency appliances that are being operated in the immediate vicinity (e.g. mobile telephones, CB radio stations, radio remote controls and microwave ovens).

Should irregular or inaccurate errors be indicated in the display ②, either remove the scale from the source of the disturbance or switch the source of the disturbance off whilst you are using the appliance.

## Cleaning

---

Clean your scale exclusively with a neutral detergent and a soft cloth.

## Troubleshooting

---

### **Warning of potential appliance damage:**

Do not attempt to dismantle the scale or to carry out repairs yourself. You could potentially damage the sensitive electronic components.

Repairs to your scale may only be carried out by qualified service personnel!

No liability will be accepted for damage resulting from improper use and/or the opening of the appliance!

If incorrect displays appear, please verify that the correct language is set. It could be that the weight is displayed in the wrong unit (kg or lb).

In cases when nothing is indicated in the display , please check to ensure the batteries are installed according to their polarities. If need be, correct the batteries insertions as appropriate.

If these measures do not correct the problem, please contact the appropriate service partner for your country.

## Disposal

---

**Do not dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance through an approved disposal centre or at your community waste facility. Observe the currently applicable regulations. In case of doubt, please contact your waste disposal centre.

### Disposal of batteries/accumulators

Used batteries/rechargeable batteries may not be disposed of in household waste. Every consumer is legally obliged to dispose of batteries at a collection site of his community/city district or at a retail store.

The purpose of this obligation is to ensure that batteries are disposed of in a non-polluting manner. Only dispose of batteries/cells when they are fully discharged.

Dispose of all packaging materials in an environmentally friendly manner.

## Warranty and Service

---

The warranty for this appliance is for 3 years from the date of purchase. The appliance has been manufactured with care and meticulously examined before delivery.

Please retain your receipt as proof of purchase. In the event of a warranty claim, please make contact by telephone with our Service Department. Only in this way can a post-free despatch for your goods be assured.

The warranty covers only claims for material and manufacturing defects, but not for transport damage, for wearing parts or for damage to fragile components, e.g. buttons or batteries. This product is for private use only and is not intended for commercial use. The warranty is void in the case of abusive and improper handling, use of force and internal tampering not carried out by our authorized service branch. Your statutory rights are not restricted in any way by this warranty.

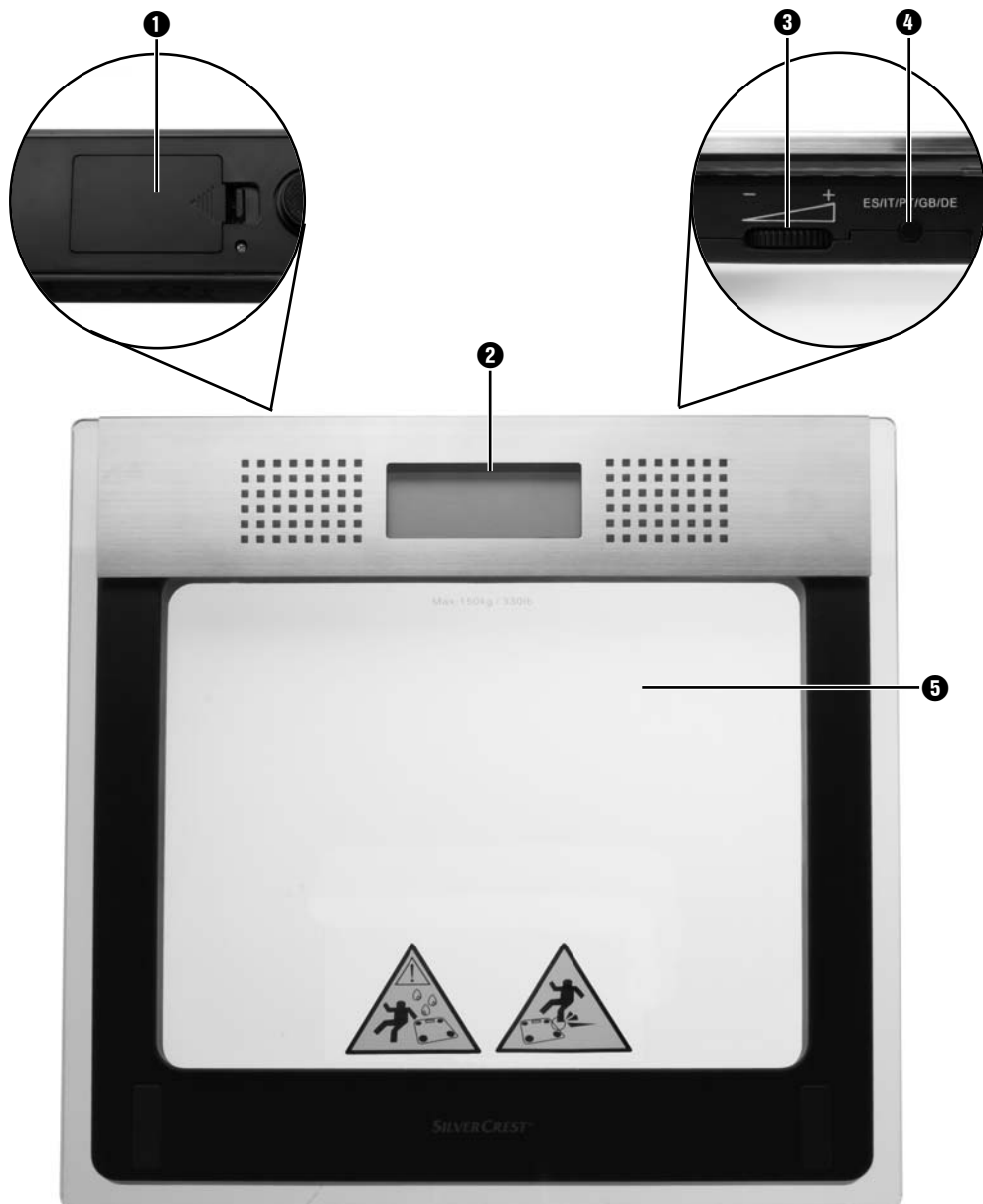
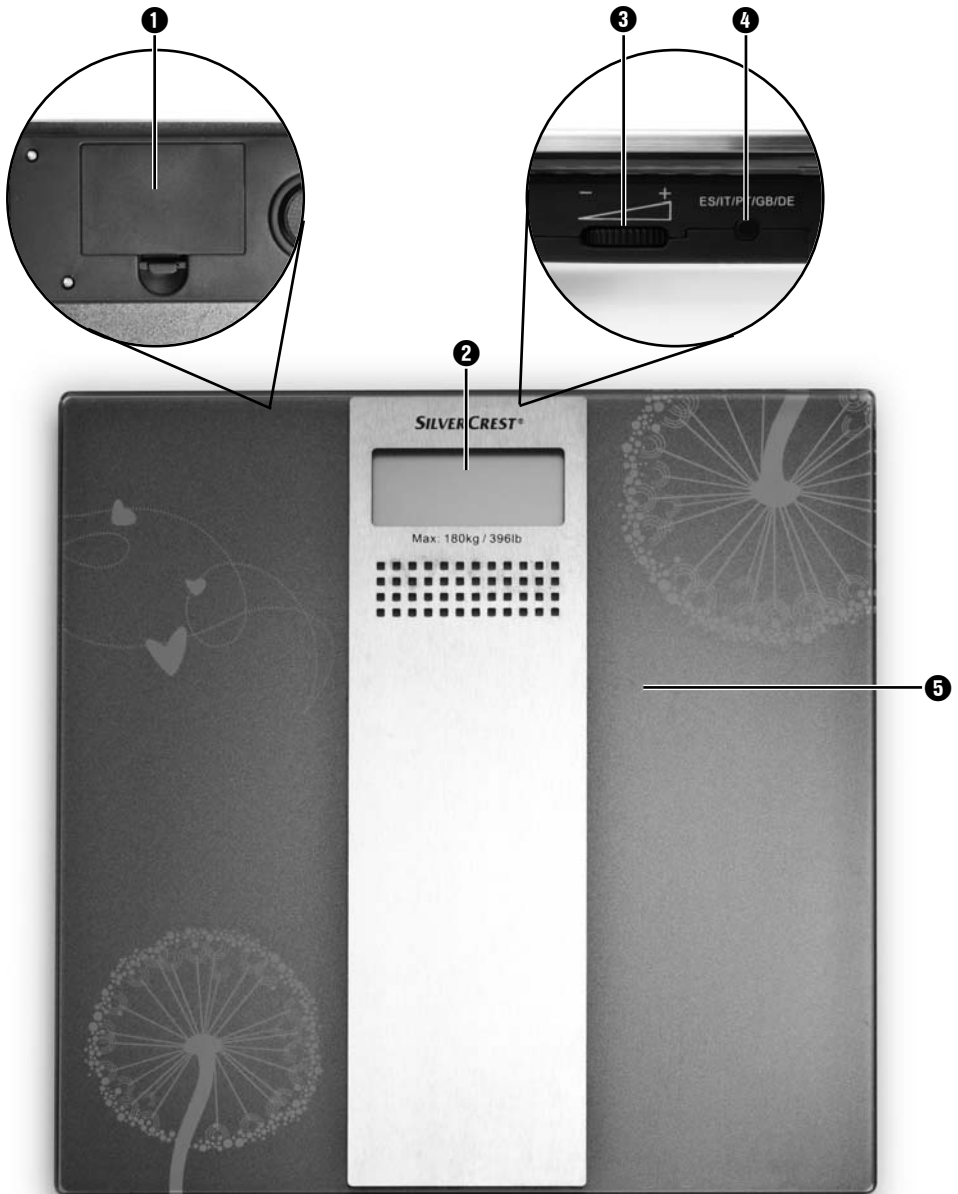
The warranty period will not be extended by repairs made under warranty. This applies also to replaced and repaired parts. Any damage and defects extant on purchase must be reported immediately after unpacking the appliance, at the latest, two days after the purchase date. Repairs made after the expiration of the warranty period are subject to payment.

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